

Doctors Friend suction units

DF-750 Suction Unit



The Doctors Friend DF-750 suction pump has a very robust ABS plastic case housing the 40W thermally protected electric motor. This model is ideally suited for doctors surgeries, small clinics and out patient departments. The DF-750 has a 1000 c.c. collection jar which is manufactured from toughened plastic to help prevent accidental breakage. An important feature is the overflow prevention valve which, along with a hydrophobic filter, gives total pump protection. All models come with a vacuum control valve. The unit weighs only 5.1kg, thus making it very portable.

**Super low maintenance
Oil-less vacuum
pump**

The DF-506, the big brother of the 702, has a larger 200w suction motor as well as a larger collection jar capacity of 2000 c.c. thus making it a perfect companion for busy nurses and doctors in small operating theatres, busy clinics and surgeries. This model has a flow rate of 40 liters per minute, double that of the smaller DF-702 unit.

Standard accessories include:

- 180cm suction tube
- Suction Catheter
- 2 filters
- 2 fuses

Optional accessories:

- Mobile upgrade as pictured below (at time or order only)

● **Optional model**



DF-506 Suction Unit



Doctors Friend suction units


DF300 Surgical and Gynae Suction

DF600 Surgical Suction

DF500 Gastric-Thoracic Drainage Pump



DF700 Mini ENT Treatment Unit

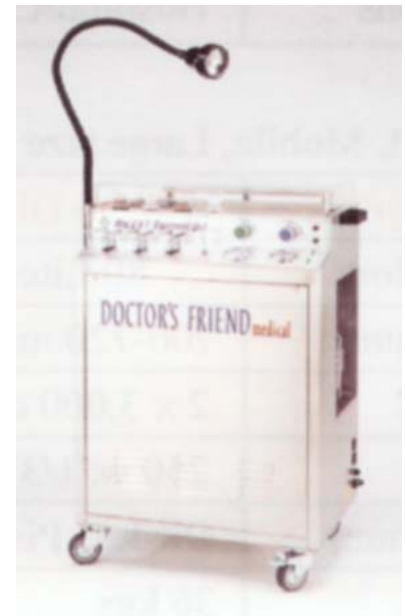

Specification	Model 300	Model 600
Vacuum Flow	80 LPM	80 LPM
Max Vacuum	700-720 mm/Hg	700-720 mm/Hg
Bottle P.C.	2 X 3000cc	2 X 3000cc
Motor	250w, 1/3 hp	250w, 1/3 hp
Vacuum Pump	Oil-less Piston	Oil-less Piston
Weight	36 Kg	36 Kg
Warranty	2 Years	2 Years
Application	Hosp/Ob/Gyn/Clin	Hosp/Sur/Clin